EN

LEVABO®

Simple solutions - Big impact

turnall® Adapters



01

Adapter for 15- or 30-degree tilt



02

Adapter for single-sided tilt or double-sided tilt/hug-function



Adapter for hug-function in all cycles

03




04

Adapter for static position or normal automatic tilt function

Easily installed
between the pump
and the hose
on the tilt system


turnall® Adapters

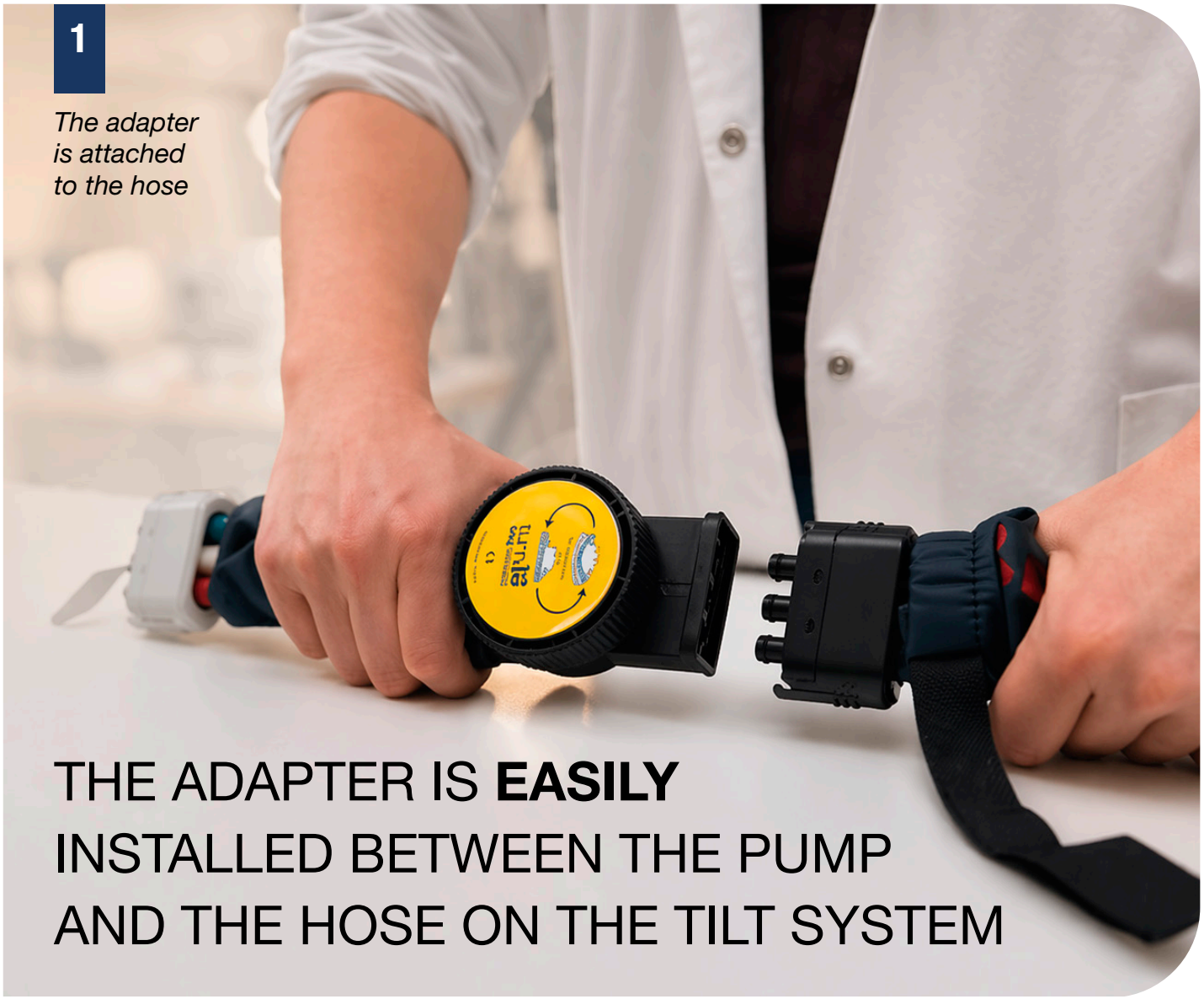

1

The adapter is attached to the hose

THE ADAPTER IS **EASILY** INSTALLED BETWEEN THE PUMP AND THE HOSE ON THE TILT SYSTEM

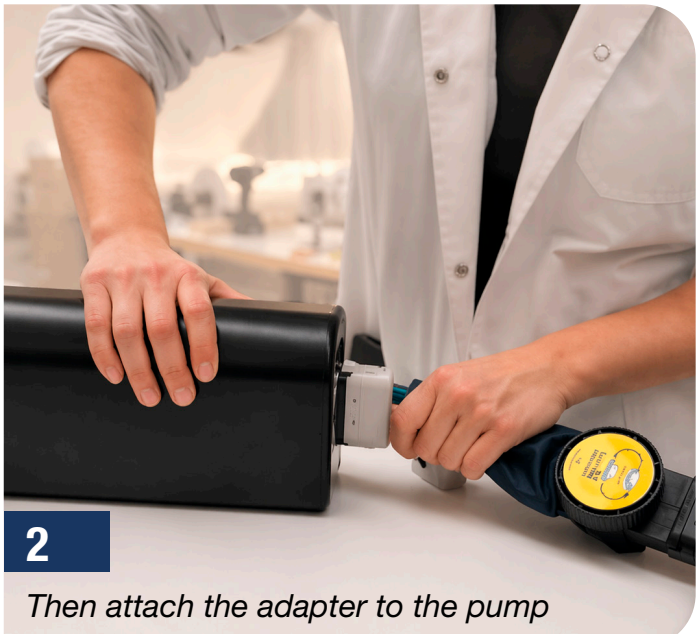

2

Then attach the adapter to the pump

You can combine both 2 and 3 adapters for additional benefits

See combination options for each adapter or in the chart on the opposite page.



2 adapters in combination



3 adapters in combination

turnall® Adapters

THE RIGHT ADAPTER PROVIDES PEACE OF MIND AND FLEXIBILITY

Choose the adapter that best suits the user's needs
and create a daily life with greater freedom, safety, and comfort

Adapter	2-in-1	Function	3-in-1	Function
01 Adapter for 15 or 30 degree tilt	🔴🟢	15 or 30 degree tilt + static position	🔴🟢⚫	15 or 30-degree tilt + static position + single-sided tilt or hug
	⚫🟢	Hug or single-sided tilt at 15 or 30 degrees		
	🟡🟢	Hug function in 15 or 30-degree position		
02 Adapter for single- or double-sided tilt/ hug function	⚫🟢	Hug or single-sided tilt at 15 or 30 degrees	🔴🟢⚫	15 or 30-degree tilt + static position + single-sided tilt or hug
	🔴⚫	Hug in static position		
03 Adapter for hug function in all cycles	🟡🟢	Hug function in 15 or 30-degree position		
	🔴🟡	Hug function in all cycles + static position		
04 Adapter for static position or normal automatic tilt function	🔴🟢	15 or 30 degree tilt + static position	🔴🟢⚫	15 or 30-degree tilt + static position + single-sided tilt or hug
	🔴⚫	Hug in static position		
	🔴🟡	Hug function in all cycles + static position		

By default, Turn All is set to a 30-degree tilt. To offer more options, we have developed four adapters. All adapters can, of course, be attached individually. In addition, they can be combined in sets of two or three. The table shows the possible combinations and their functions. *Note that in combinations involving the red adapter, it must **always** be placed closest to the pump.* The color combinations must be followed and attached in the order shown, from the pump to the hose.

turnall® Adapters

01 | Adapter for 15 or 30-degree tilt


Item No. 76167



The green adapter is used when a 15-degree tilt is desired instead of the preset 30 degrees that Turn All normally operates at.

Who uses it?

The adapter may be relevant for individuals who are sensitive to changes in their daily routine or skeptical of repositioning and positioning. In such cases, a 15-degree tilt can be a gentler solution. In practice, this means a 15-degree tilt on one side and approximately a 5-degree counter-tilt on the opposite side. If you wish to return to the standard setting, simply remove the adapter or rotate it to 30 degrees.

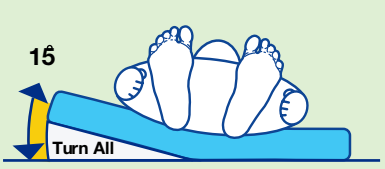

Position A

15-degree tilt (in all cycles)



Position B


30-degree tilt (in all cycles)





02 | Adapter for single-sided tilt or double-sided tilt/cradle function

Item No. 76210

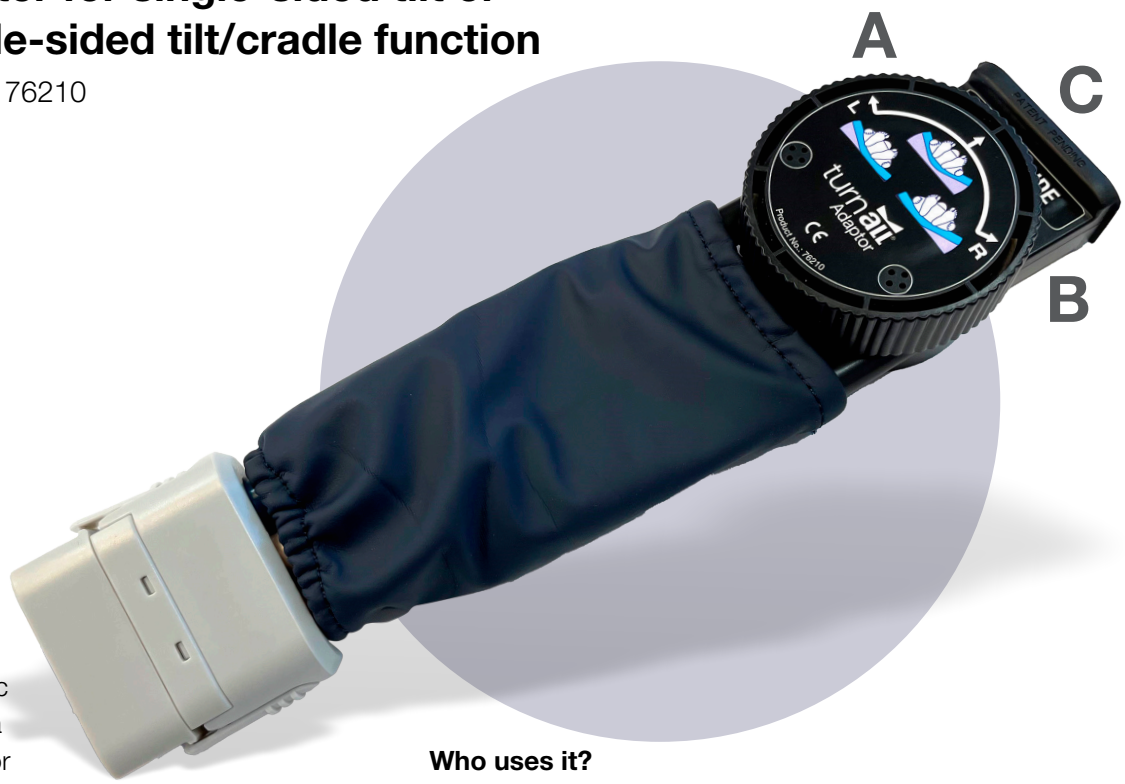

The black adapter
is used when tilting is desired in only one specific direction. For example, if a patient does not want to or is unable to lie on one side, the adapter can be adjusted to tilt exclusively in the desired direction. *(Turn the round wheel.)* The adapter can also be set to a “hug” function, where the patient is not simply “tilted” from side to side, but instead lies securely in an enveloping position. This can be beneficial for patients who need increased comfort. *(Turn to the symbol for the “hug” function.)*

Who uses it?
For single-sided tilting, the adapter may be suitable for individuals with, for example, pressure ulcers or fractures on one side of the body, or those with a preferred lying position. The hugging function may be suitable for individuals with scoliosis, as well as those who need increased security in their daily lives. In addition, it can be beneficial for residents who need enhanced body awareness, as well as those with mental health conditions.




Position A
Single-sided tilt (left)




Position C
Hugging function / double-sided tilt



Position B
Single-sided tilt (right)

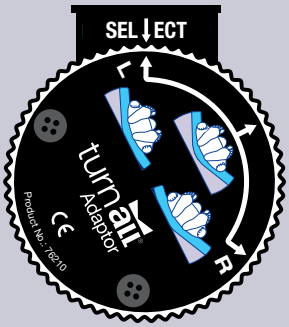

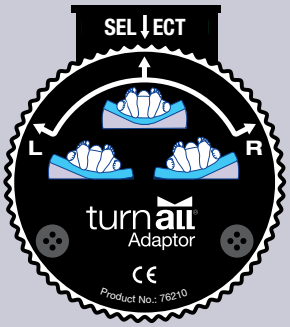



03 | Adapter for hug position in all cycles

Item No. 76168




The yellow adapter is used when you want to lie in the hug position for the time interval set on the pump (30, 60, or 90 minutes). If, for example, 30 minutes is selected, the patient first lies in the hug position for 30 minutes and then on their back for 30 minutes. After that, the system automatically alternates between the “hug” position and the supine position at the selected interval. Positioning in the “hug” position is beneficial because it promotes a sense of security, calm, and physical well-being through the release of the hormone oxytocin. It is an effective method for reducing stress, alleviating pain, and improving sleep quality.

Who uses it?

The adapter may be suitable for individuals who have had a blood clot, are paralyzed following an accident, are post-operative, or are temporarily immobile due to medication or a coma.

Hug function

Double-sided tilt as follows:
30, 60, or 90 minutes.

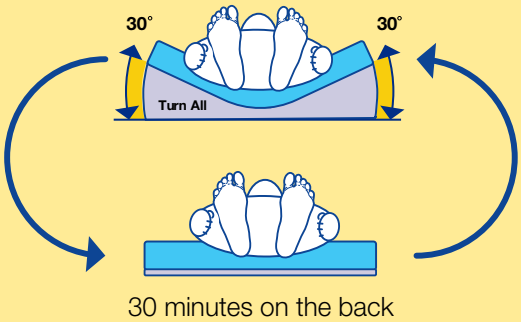


30 minutes on the back



04 | Adapter for static position or normal automatic tilt function


Item No. 76175

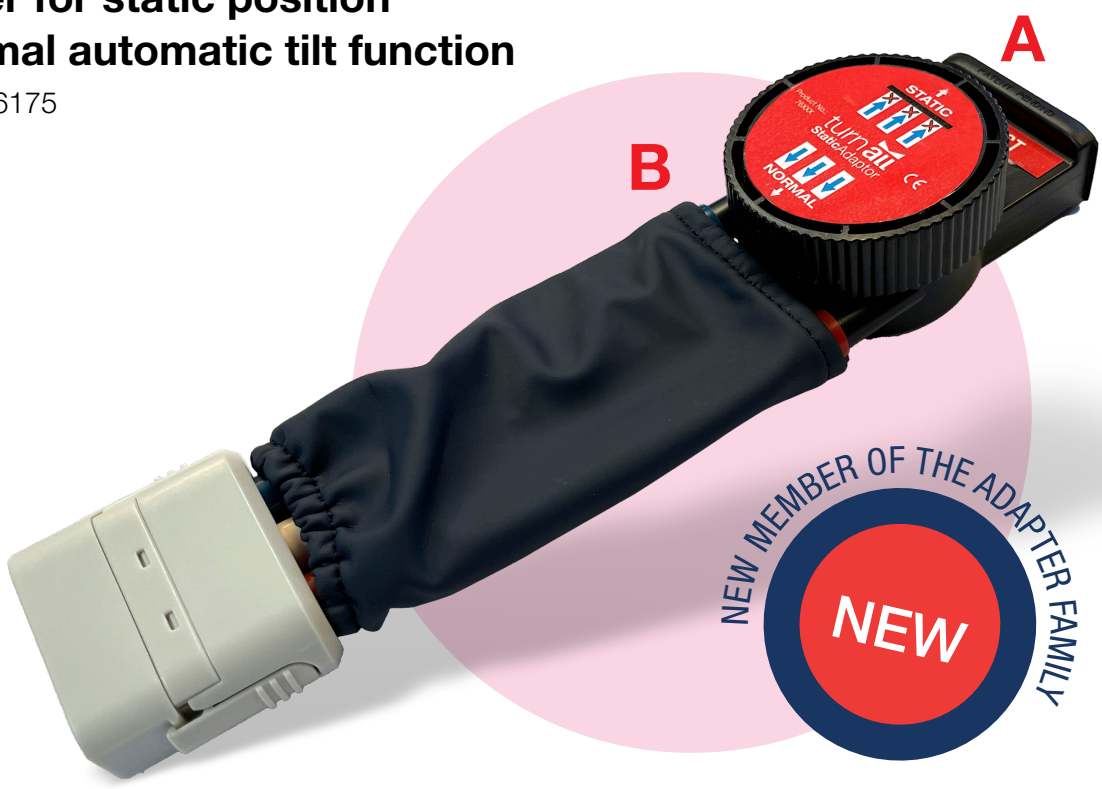

The red adapter is used when a static tilt to one side is desired without an automatic tilt function. It provides stable and controlled positioning, allowing the patient to remain in that position until a new position is selected.


Who uses it?

The adapter is particularly relevant for patients who require calm and stable positioning. This may include, for example, patients with numerous tubes, drains, or catheters, where movement can increase the risk of dislodgement. It can also be used for patients with respiratory challenges, where a stable position supports breathing. In addition, it may be relevant for patients in pain or those who are restless, who may react negatively to repeated movements, as well as for immobile or frail patients, where static positioning is important to prevent complications such as pressure ulcers and discomfort.

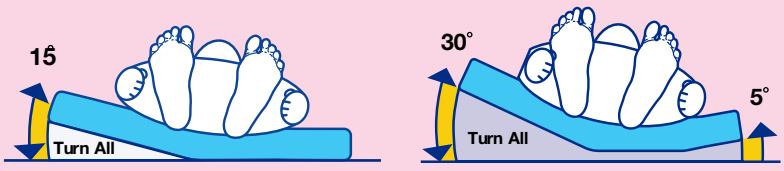

Position A
Static position

Position B
Normal tilt function

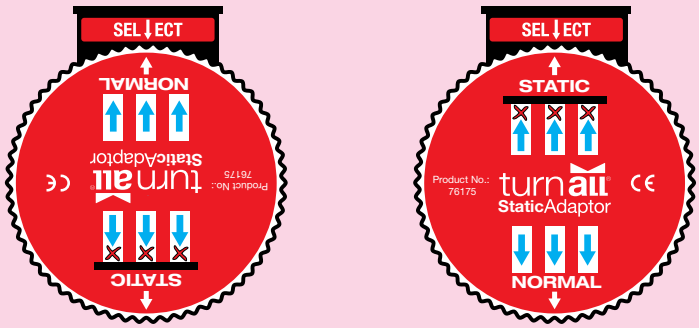

Welfare technology
that helps **you**
help **others**












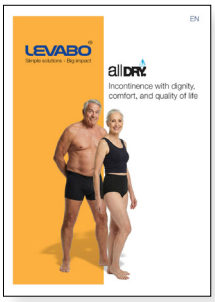



Scan the QR code
to access our catalogues for All Dry,
All Up, Heel Up, Seat All, Turn All,
Adapters og MultiPatient

Levabo is certified according to ISO-standard 13 485

UK CA CE All Levabo products are CE marked
and fully comply with MDR 2017 / 745

ISO 13485
BUREAU VERITAS
Certification

LEVABO®
Simple solutions - Big impact

HOSPITAL HOME CARE NURSING HOME

Levabo ApS
Sverigesvej 20A
8660 Skanderborg
Tel.: +45 31 70 40 25
Mail: info@levabo.dk
www.levabo.com

Turn All Adapter Total EN 28.10.2026 V4